



User Guide




Chapters

1. Overview of features
2. Your settings
3. Scanning
4. Switch access
5. Dictionaries
6. Voices
7. Application Themes
8. Keyboard Layout
9. Categories
10. Adding A Category or Phrase
11. Editing Phrases
12. History
13. Backing Up
14. Email
15. Facebook & Twitter
16. SMS
17. Emoticons
18. Hand writing input
19. Comparison Table – Functionality by iDevice using various access methods

CHAPTER 1

OVERVIEW OF FEATURES

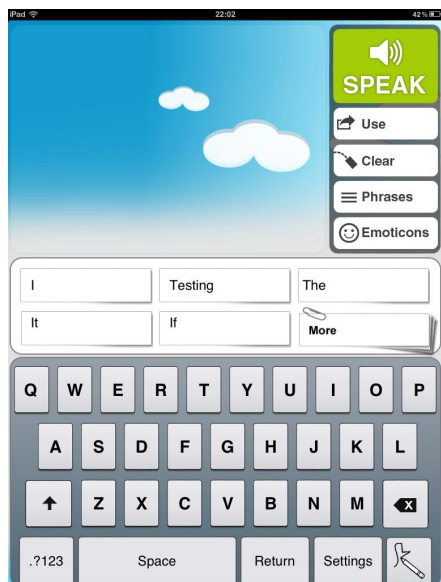

Voice Output

The user is able to speak out any message that is input into the message window or from the Categories pages or History list. The user is able to choose from a range of Loquendo™ voices and emoticons.

Access Methods

There are three available access methods for using Predictable. These are direct access, where the user types using the on screen keyboard (choice of QWERTY, ABC, high frequency layouts); “Touch Anywhere”, where the user makes use of the scanning function and makes their selection by touching anywhere on the screen to select the section highlighted on the scan. The third method is access via a switch(es). **This requires a Bluetooth switch interface** (Therapy Box Switch Box, switch interface from RJ Cooper, Zyteq, Spectronics or Compusult). The user’s switch(es) is plugged into this box. This requires syncing the Switch Box to your iPad, iPhone, iPod Touch. An assistant would then enter their preferences for scan pattern to best suit their needs. If using two switches, the user can select “User Scan” mode to allow for a scan + enter method. These three modes can be used independently only; and not in conjunction. It is simple to switch between modes.

Prediction

Predictable currently features two dictionaries, a UK and US based dictionary. The

prediction occurs based on letters entered and also on next word prediction. The display brings up 5 predictions on screen as well as a tab for more, which then brings up a further set of 5 predictions. The prediction is influenced by user pattern and will “learn” new words entered and the context in which these words are used in. The prediction engine also manages high frequency spelling errors (e.g. predicting “the” when “hte” is entered).

Customisation

Users can customise the look and feel of the display by choosing from two themes. These can be changed as the user likes. The user can also choose from a range of voices. These voices are from Loquendo™. Users can choose from British, American and Australian voices. Users can switch between voices if they wish.

Customisation of the keyboard layout is also possible.

Categories (formerly called My Phrases)

This page organises a users phrases into quick access categories. The app comes with some suggested categories, but can be easily edited to suit the user. Each category includes a list of words or phrases which can also be deleted, added or edited. Each item can have an image attached to it. There is a simple on screen editor which allows users or their carers to complete this quickly.

Adding A Phrase

Users can add a phrase from items typed in the message window or History page or they can add a phrase directly into the Categories section.

History

This section provides quick access to the most recently entered words or phrases placed (and spoken or published) in the message window. This can be accessed to repeat something recently said, or can be used by carers, teachers or therapists to help gauge how someone is using the application. This feature can also be turned off for privacy reasons.

Backing Up

It is possible to back up content of Categories and History to a computer to save. These can also be returned to the application after changes are made.

Sending an email

It is possible to write the content of an email and integrate it with the email program

you have running on your iPad, iPhone or iPod Touch when using wifi, or your 3G network.

Sending an SMS

This is only possible if you are using the application on an iPhone while using Direct Touch access method. You are able to integrate the content of your message window into an SMS.

Updating your Facebook status or sending a tweet

You are able to update your status on Facebook or send a tweet on Twitter using the text prepared in your message window. Please note this is only available when using it in direct typing access.

CHAPTER 2

YOUR SETTINGS

It is important to set up Predictable to suit the user who will be using Predictable.

The Settings are available outside of Predictable. To access these Settings:

1. Go to Settings on your iPad/iPhone/iPad Touch




2. Scroll down to Predictable



Settings – Choose Access Method

Choose one of the following:

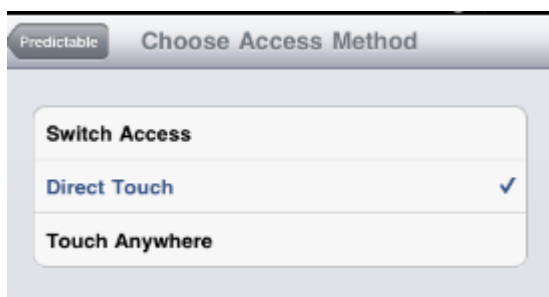

Direct Access - for using your finger or stylus to type messages on the on screen keyboard or note pad. Not for use with scanning or switches.

Switch Access – for use with one/two switches. Remember switches cannot be plugged in directly to an iPad/iPod/iPhone and require an intermediary. This application works with the Therapy Box Switch Box. You will need to pair your Switch Box with your iPad/iPod/iPhone.

Pairing your Bluetooth Switch Box with your iPad/iPod/iPhone:


1. Go to Settings on your iPad/iPhone/iPad Touch



2. Scroll to “General”



3. Turn Bluetooth “On”

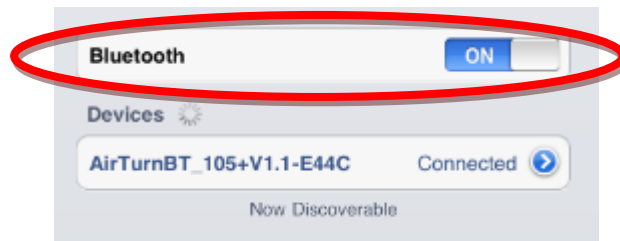


4. Plug in your switch to the switch box

5. Turn on your Switch Box by pressing the red power button and holding for 1 second.

The LED will blink 3 times and attempt to pair with your iPad/iPod/iPhone. The LED will blink orange rapidly while it is connecting. Once connected, it will blink rapidly green and red alternately. Then release.

6. Choose your Device

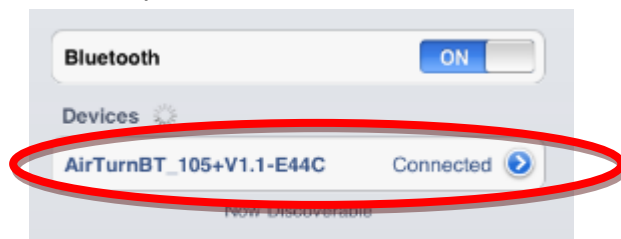

7. Wait for it to find the Device and display “Connected”

If pairing is unsuccessful, delete the device displayed in the host. Reset and try pairing again. Note, initial discovery may take up to 30 seconds.


To reset pairing (e.g. to set up with a different iPad/iPhone/iPod), press and hold the red button for 8 seconds until it blinks rapidly green and red alternately. Then release.

Please download the switch user guide from www.tboxapps.com for further detailed guide on switch integration and usage.

Touch Anywhere - In choosing this, the scanning mode will default to auto-scan mode.

Settings – Scanning Type

There are two scanning modes you can choose.

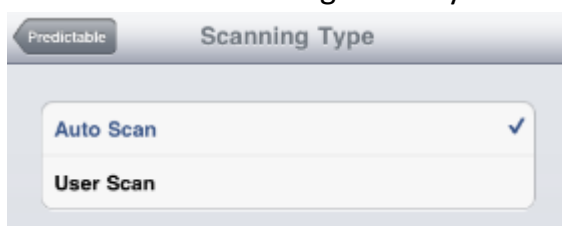

These are:

1. Auto Scan – This scans each section one-by-one. You can adjust the scan rate and the number of scan cycles.
2. User Scan – this is sometimes known as “tap to advance” and requires two switches to operate in this mode. In this mode, the first switch moves the scanning highlight between sections, and the second switch is used to select the required section.

Please note Auto Scan can only be used with either of 1) one switch or 2) Touch Anywhere mode. User Scan can only be used with two switches.

Settings – Scanning Rate



Choose the speed that you would like the scanning to occur. There are five settings (very fast, fast, normal, slow, and very slow).

Settings – Stop Scanning



Choose the number of times you would like the scan pattern to continue. This may be useful if the user requires extra time to think /find their selection. You can choose between one to five scan cycles. You can also choose unlimited scans.

Settings – Force Delay

Enabling the Force Delay option will mean that a set delay period is activated before the next tap on the screen is registered. The options available are:

Off

1 Second

3 Seconds

5 Seconds

Settings – Mail Settings

To use the email settings within Predictable you need to set up your email account in the settings. Choose from the following email hosts:

Yahoo

Hotmail

Google

Mobile Me

Others

You will then need to enter your email address and password to be saved in the settings. If you have an email host other than the four listed above, you would need to enter the host smtp address.

You can also enter a default To email address if there is one email address that the person using Predictable would like to primarily send emails to. Within the app, they can choose this default address or search through their other contacts.

Settings – Keyboard Type


There are three available keyboard layouts:

QWERTY

ABCD (alphabetical order arrangement)

High Frequency layout

Settings – Dictionary

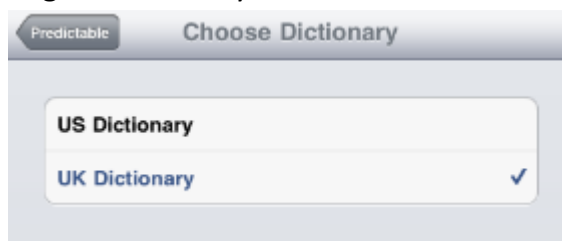

Choose either the UK or US generated dictionaries. These will impact on the prediction options and also the spelling used by the prediction engine (e.g. colour vs. color).

Settings – Speak or Publish

The default setting for Predictable is for Speak. Publish is another option if you do not wish to have a message spoken. It will allow you to type your message and then tap Publish to invert your message for a communication partner sitting opposite to read.

Settings – Save History

This setting can be turned on/off. If you choose “off” the History page will not be populated with the messages typed in after choosing this setting.

Settings – Access Emoticons

This setting can be turned on/off. If you choose “off” the emoticon function will not be available when using Predictable.

Settings – Speak When Highlighted

This setting can be turned on/off. When enabled, this will provide a spoken

representation of the button within Predictable including keyboard keys, prediction options, function keys etc. There are some items on the keyboard that it does not currently speak aloud which are some punctuation keys.

CHAPTER 3

SCANNING

Scanning is only required if you wish to use a switch(es) or have chosen the Touch Anywhere mode. Refer to Settings information to make sure you have set up as per your requirements.

The scanning speed can be adjusted within Settings. There are five scanning rates:

Very Slow

Slow



Normal

Fast


Very Fast

You can also adjust how many scanning cycles should take place if the switch is not activated.


The scanning works by highlighting each of the key working sections one-by-one. When the highlight is on the area you need, you need to press your switch (if using a switch) or touch the screen (if using Touch Anywhere). The scan will then highlight the items within that area, for example, the keyboard letters.

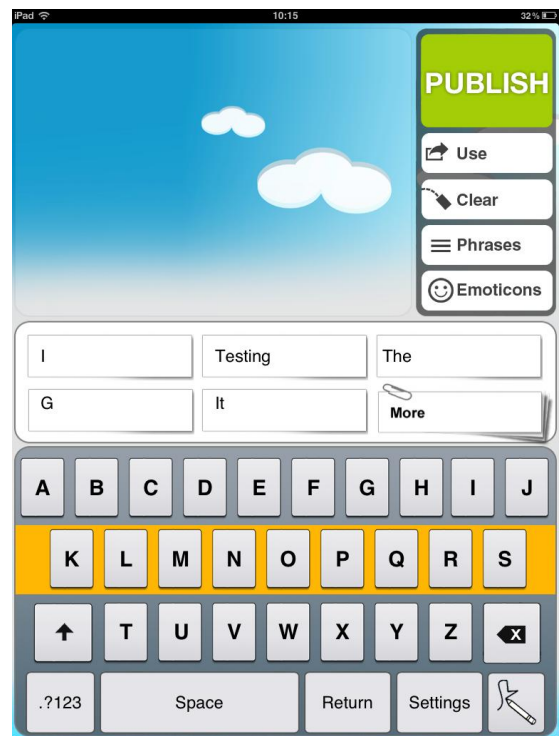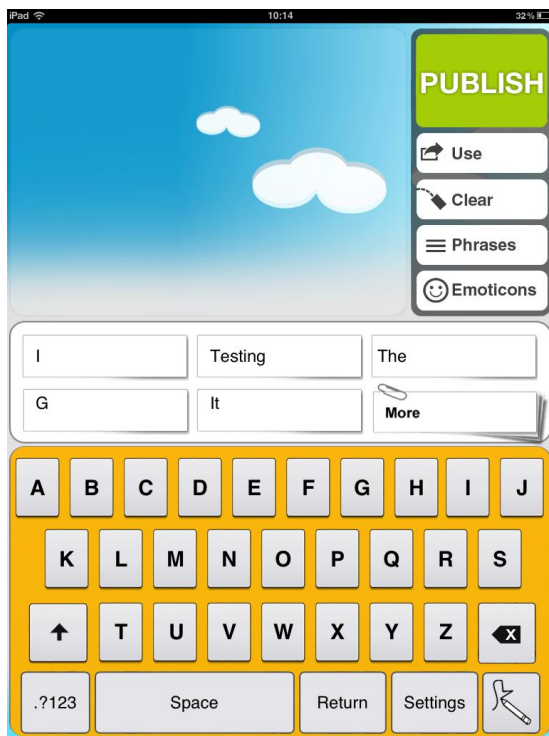


When scanning the keyboard, the scan will cover rows and then the individual columns.

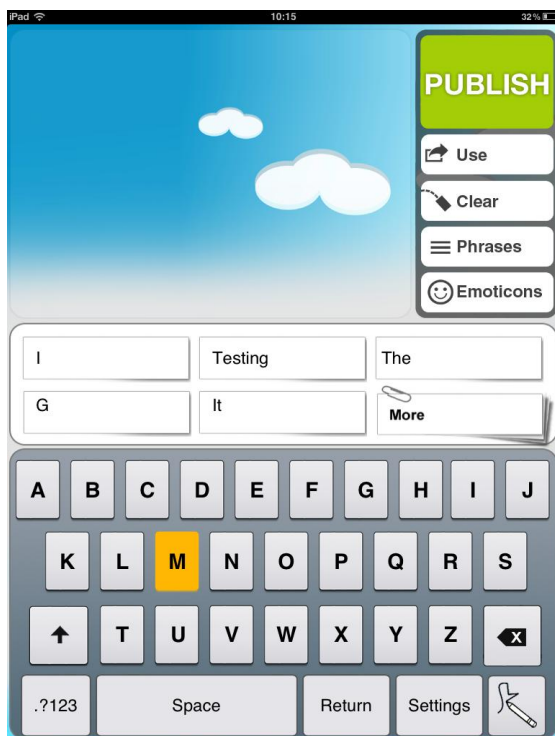

When scanning the message window, and the Speak (or Publish), Use, Clear, Phrases and Emoticons buttons, it will scan the message window first and then go through the five buttons to the right.

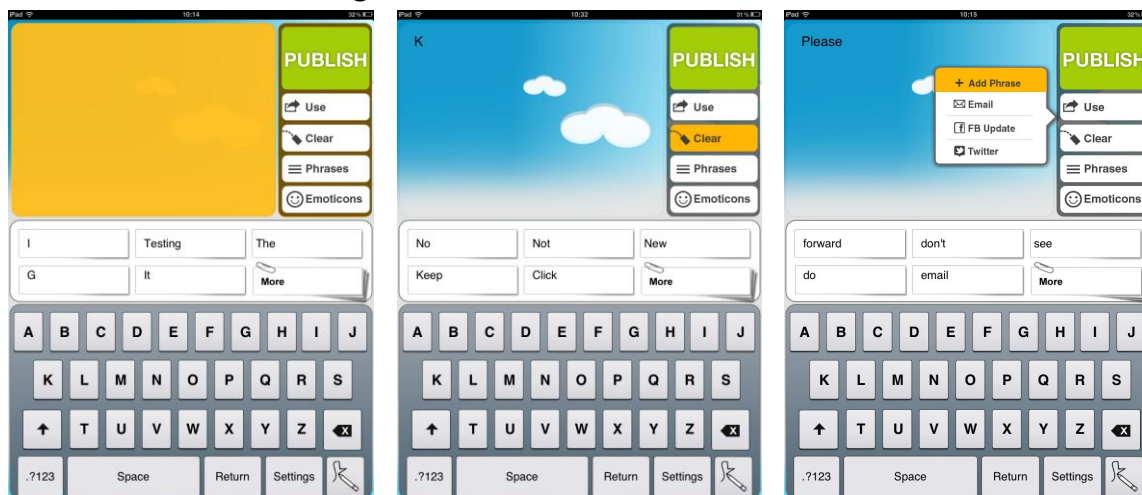

When scanning the prediction options, it will scan the rows and then the three items within each row. You may choose to have auditory feedback of each of the predicted items scanned. This feature is turned on/off in the Settings.

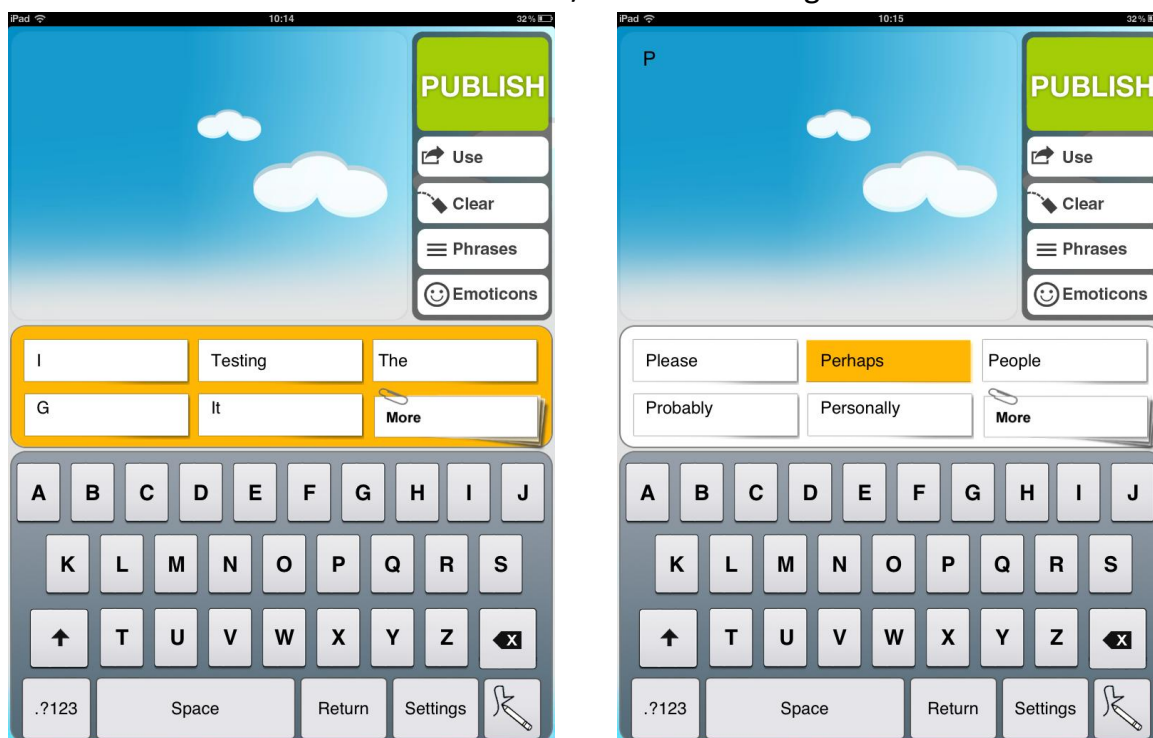

When scanning in Categories, the scanning will again highlight each of the working sections including the back and forward arrows to move between pages if there are multiple pages of categories.



When scanning the Categories folders, it will scan each row, and then when a row is chosen, scan each folder across that row from left to right. If a particular row is not chosen it will then scan the back and forward arrows.

When scanning the phrases within the Categories, it will scan each phrase and then the back and forward arrows to access further pages.



When scanning the History page, it will scan each working section. When scanning the History word /phrase list, it will scan row by row. Once a row is chosen you can then scan and select the phrase to speak/publish or scan and select to add the phrase to the Categories pages or to delete the phrase. It is not possible to use the scanning features when adding photos. You may require someone to assist you in setting this up.

CHAPTER 4

SWITCH ACCESS


To use a switch with Predictable, a Switch Box is required. The Switch Box is a Bluetooth transceiver and needs to be paired with your iPad, iPod or iPhone. Please check www.tboxapps.com for a list of retailers which supply the switch box or an equivalent.

Pairing your Bluetooth Switch Box with your iPad/iPod/iPhone:


1. Go to Settings on your iPad/iPhone/iPad Touch



2. Scroll to "General"

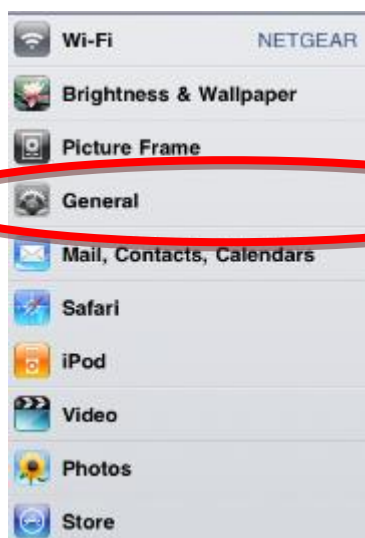


3. Turn Bluetooth "On"

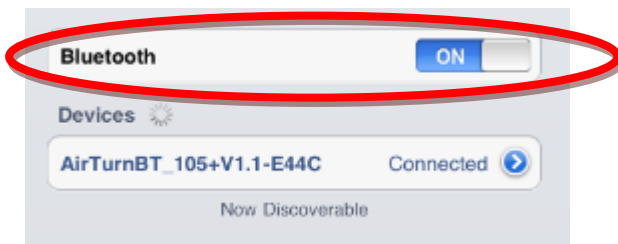

4. Plug in your switch to the switch box
5. Turn on your Switch Box by pressing the red power button and holding for 1 second.

The LED will blink 3 times and attempt to pair with your iPad/iPod/iPhone. The LED will blink orange rapidly while it is connecting. Once connected, it will blink rapidly green and red alternately. Then release.

6. Choose your Device

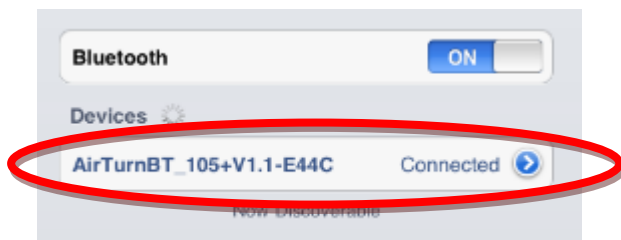


7. Wait for it to find the Device and display “Connected”


If pairing is unsuccessful, delete the device displayed in the host. Reset and try pairing again. Note, initial discovery may take up to 30 seconds.

To reset pairing (e.g. to set up with a different iPad/iPhone/iPod), press and hold the red button for 8 seconds until it blinks rapidly green and red alternately. Then release.


When using switches, you need to check you have chosen the right Access Method.

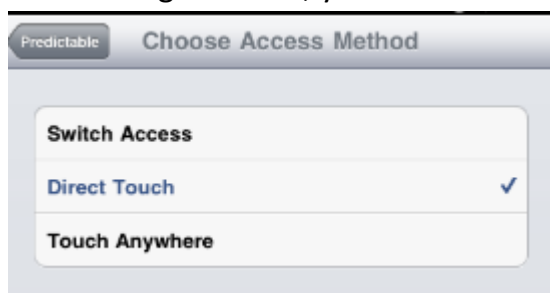

For single switch use, choose “Auto Scan” and then choose your desired scanning rate; and the number of cycles for scanning.

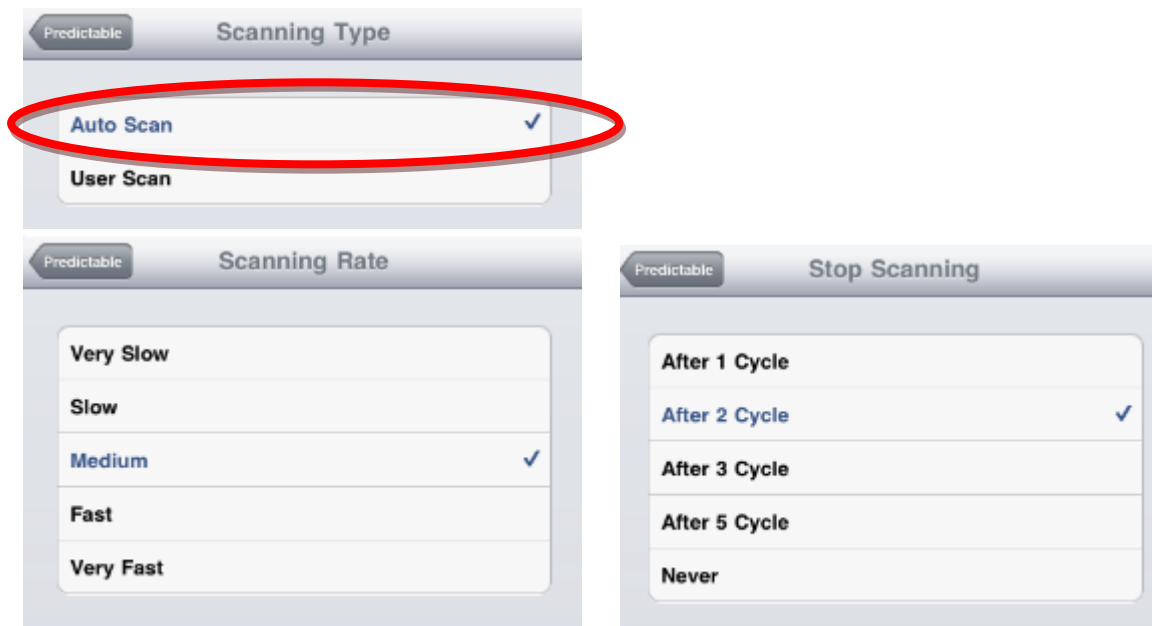


For two switch use, choose “User Scan”.

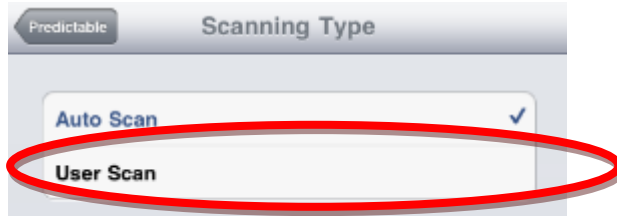

The scanning rate and number of cycles is not applicable in this situation.

For Touch Anywhere, choose “Auto Scan” and then choose your desired scanning rate; and the number of cycles for scanning.

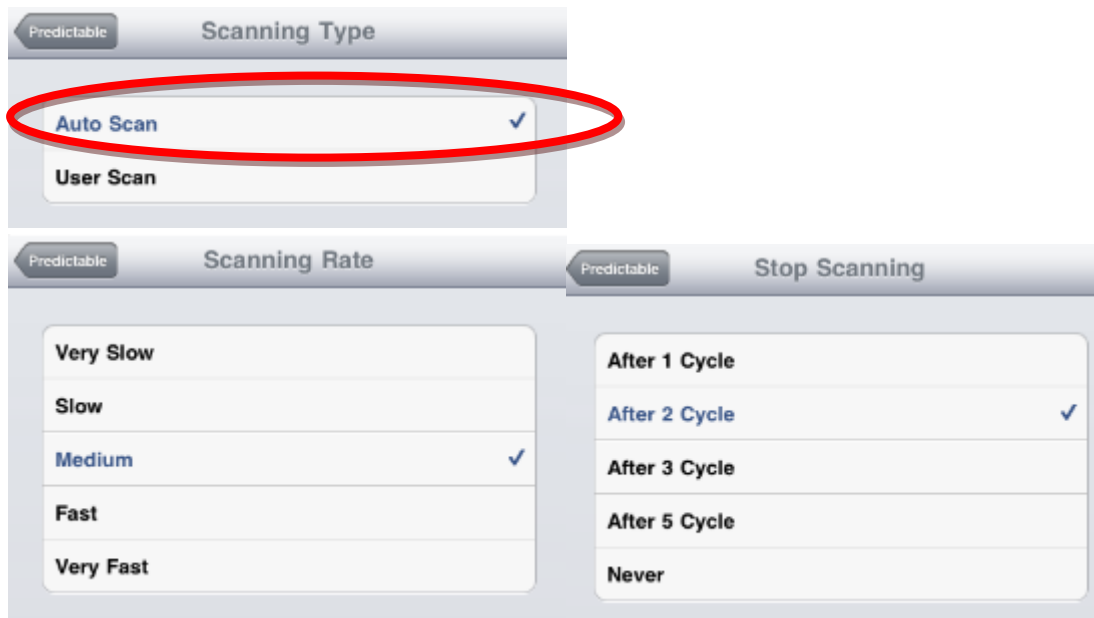

CHAPTER 5

DICTIONARIES



There are two inbuilt dictionaries, UK and US. The differences between them reflect variations in word use and selection. The spelling may not always reflect UK/US spelling rules.



The dictionary will begin to learn new words and also patterns of use. A new word or pattern of use will be “learned” within approximately two entries of either the new word or use of the same sentence/context to learn a new pattern.

Sometimes new or existing words require pronunciation changes. This can be easily done by going to the User Dictionary in the Settings within Predictable. Here you can add the word as it is spelled e.g. “iPad” and then the way it is to be spoken “eye pad”. To add or edit a pronunciation:

1. Go to Settings in Predictable
2. Go to User Dictionary
3. Tap on Add
4. Type in the way the word is spelled
5. Type in the way the word should be said.
6. Tap “Speak” to check that it is right.

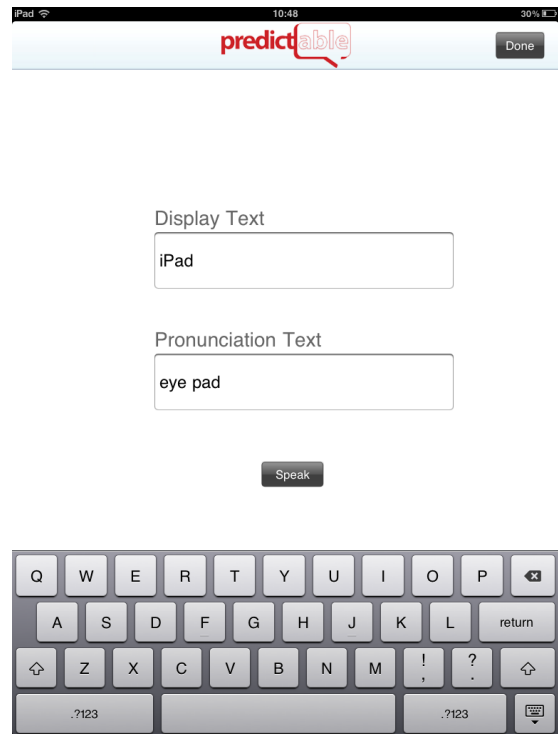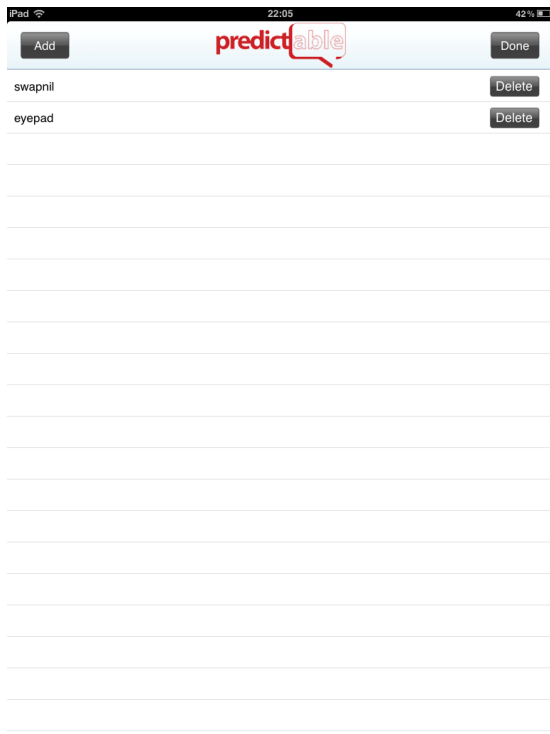

You may also wish to use this feature to add abbreviation expansions or short cuts.

For example you might type in

Display Text: btw

Pronunciation Text: by the way

Or

Display Text; QQQ

Pronunciation Text: Hi. My name is Rebecca and I live in London. Nice to meet you.

CHAPTER 6

VOICES

There are nine voices available to use. These are included with the application and do not require further downloads/purchase.

The following voices are available:

Australian

Female – Grace

Male – Alan

American

Female – Allison

Female – Susan

Male – Dave

Male - Steven


British

Female – Kate

Female – Elizabeth

Male – Simon

From within Predictable, go to the Settings icon at the bottom of the page.

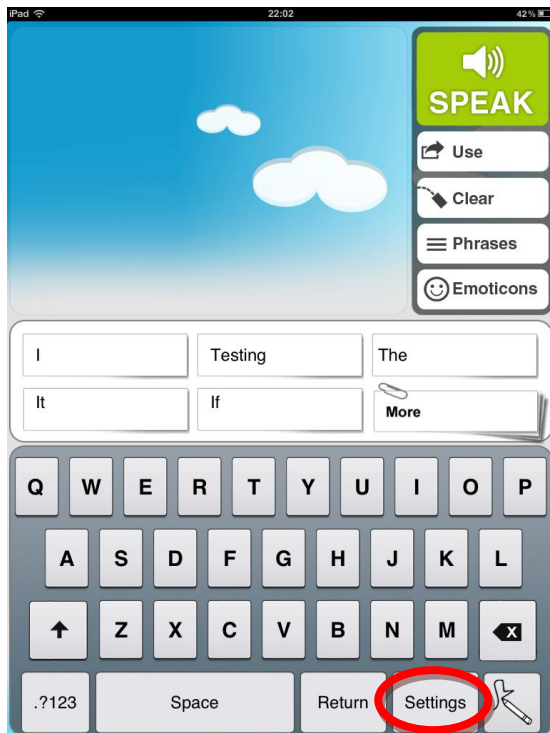

Choose “Voices”

Select your preferred voice.



CHAPTER 7

APPLICATION THEMES

From within Predictable, go to the Settings icon at the bottom of the page. Choose “Application Themes”. There are two themes currently available. Select a theme. This will act as wallpaper for Predictable.

CHAPTER 8


KEYBOARD LAYOUT

There are three different keyboard options:

- 1) QWERTY
- 2) Alphabetical Order
- 3) High Frequency Layout

To change the layout:

1. Go to Settings on your iPad/iPhone/iPad Touch
2. Scroll down to Predictable

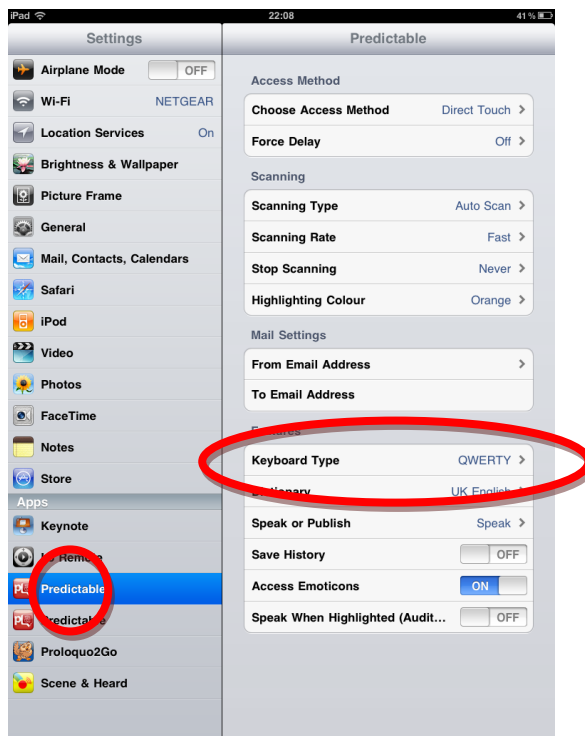


3. Choose “Keyboard Type”. Select your preferred option.

CHAPTER 9

CATEGORIES – FORMERLY CALLED MY PHRASES

The Categories page is designed as a way to store words or phrases that you can access easily and quickly. This might be words or phrases you use frequently or would like easy access to. Predictable comes ready with some example categories and phrases, but the idea is to customise the categories and phrases to suit the person using Predictable. Both the categories and phrases can be allocated a photo/symbol/image if desired.

If you have enabled “Speak when highlighted” then the phrases will be said as you tap on each phrase. If you have not then the phrase will be added to the message window to speak from there.

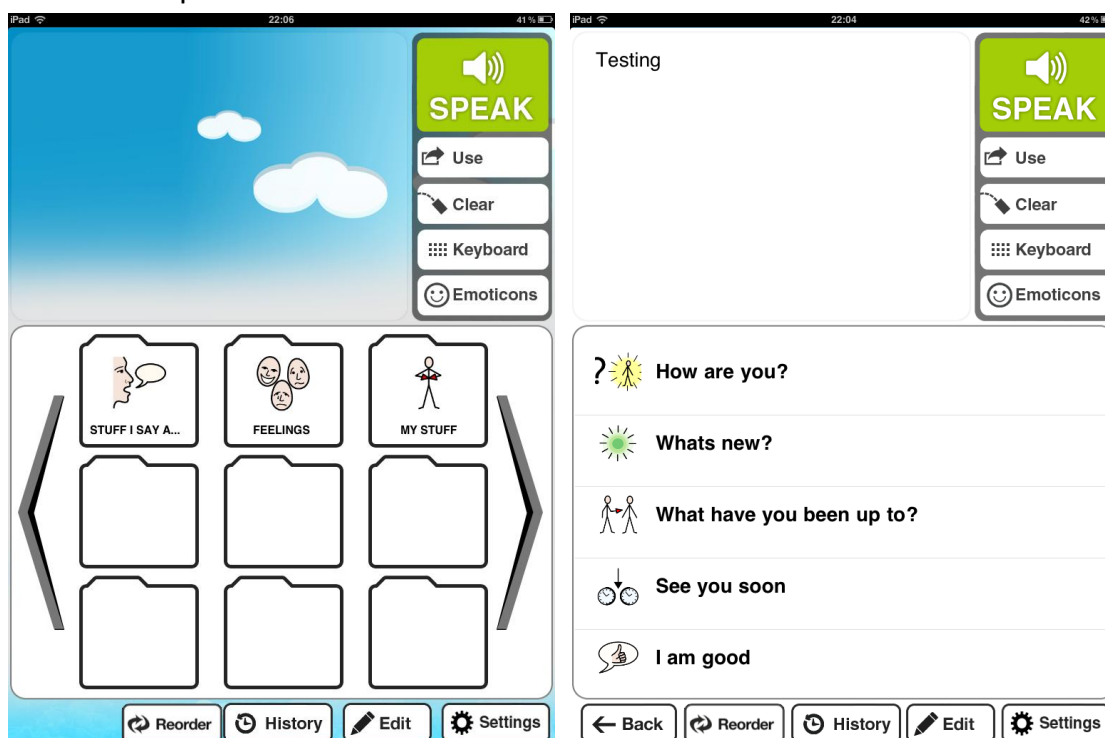

CHAPTER 10


ADDING A CATEGORY OR PHRASE

By selecting the Use and then Add Phrase button to the right of the message window, it will prompt you to add the phrase to the Categories section by asking you to select a category to add it to.



Add Phrase can be accessed in direct access, touch anywhere and scanning modes. In scanning mode, wait for the scan to include the top line (message window; speak; use; clear; phrases; emoticons); select this space, and then wait for the scan to highlight “Use” before then scanning through the Use menu to Add Phrase.

When you choose “Add Phrase” you will be prompted to choose which category you want to add it to (“Select Category”).

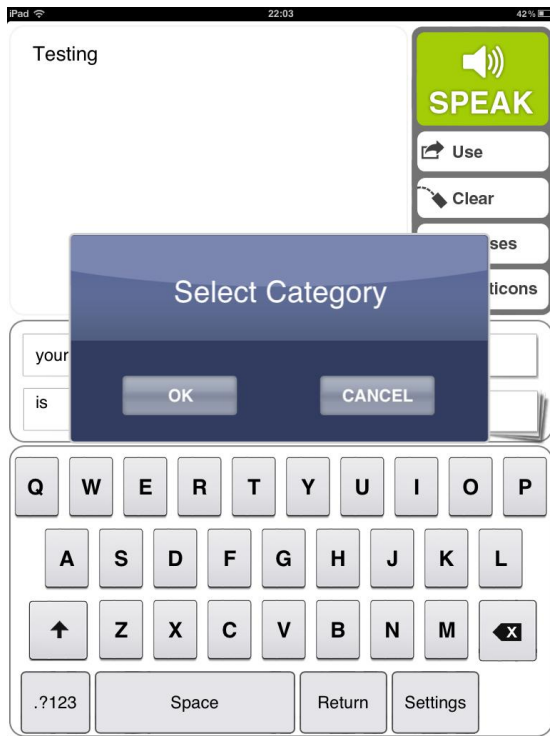


The list of categories in My Phrases will then pop up and you can choose the category you want to add it to.



You will be prompted to “Add Photo”. When you select this, your iPad/iPhone/iPod Touch photo album will appear and you are able to choose the photo you would like to allocate to the phrase. This can only be accessed in direct access mode.


You can edit the names or photos on the Category pages by tapping “Edit” and then the category you would like to edit.

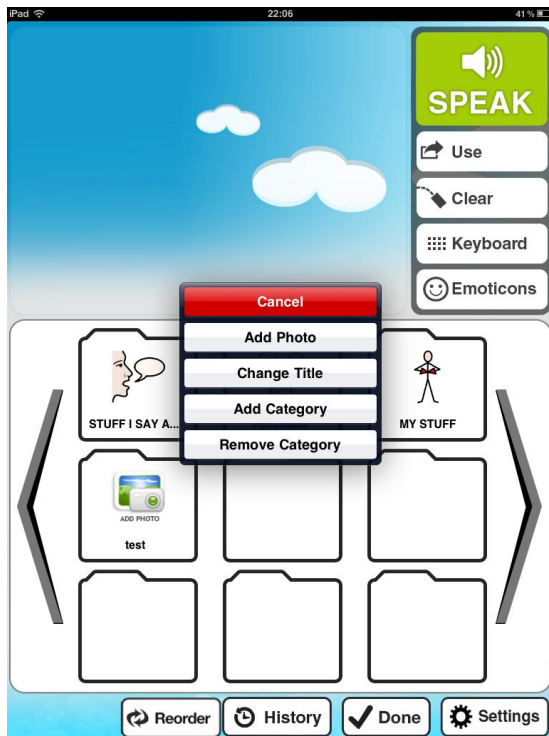

Then you are prompted with four options (cancel; add photo; change title; add category; remove category”).

To add a new phrase within the category, tap on Edit once you are in the category you want to add the new phrase to. Tap anywhere on the phrase list and you will be prompted with the four choices of Add Photo, Add Phrase, Edit Phrase, and Remove Phrase.


Choose Add Phrase – a new window will pop up with a field to enter the new phrase. Once complete, tap on Done. You can then further edit the phrase to add photos etc.

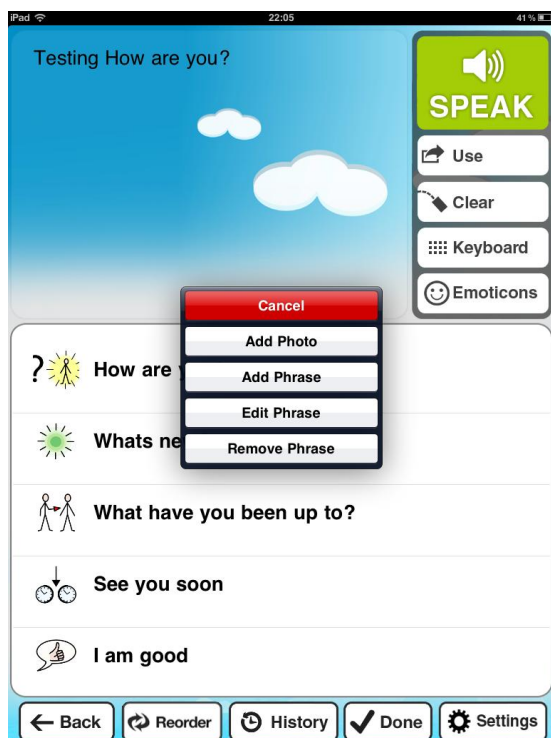

CHAPTER 11



EDITING PHRASES

When you are in the phrase list within the category:

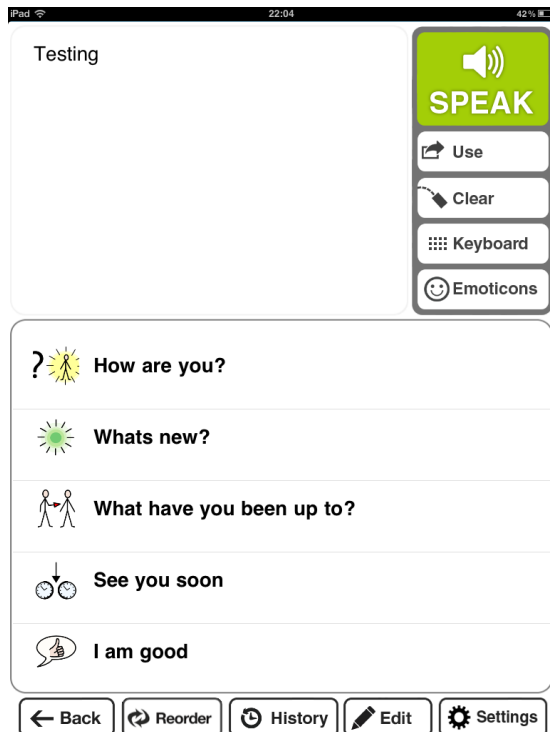

1. go into edit mode by tapping on Edit at the bottom of the screen.
2. tap any phrase or any free space in the list. It will give you five options (cancel; add photo; add phrase; edit phrase; remove phrase).

If you choose add photo, your iPad/iPhone/iPod Touch photo album will appear and you are able to choose the photo you would like to allocate to the phrase. If you choose edit phrase, you can re-type the phrase using the message window.

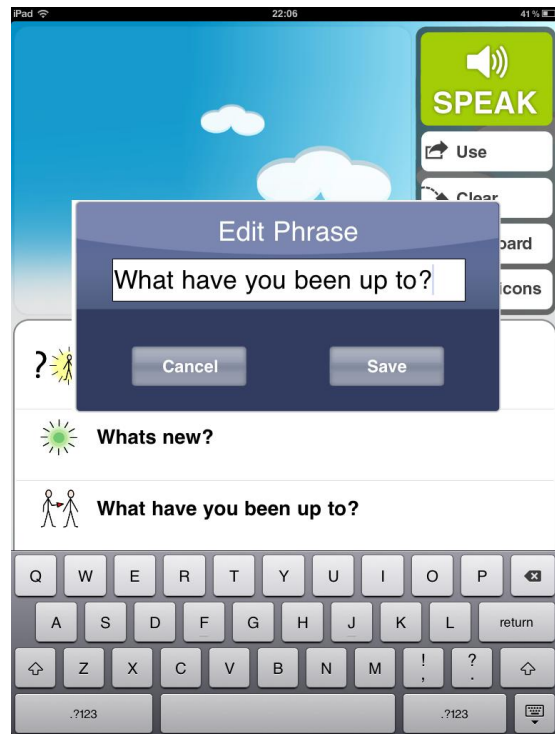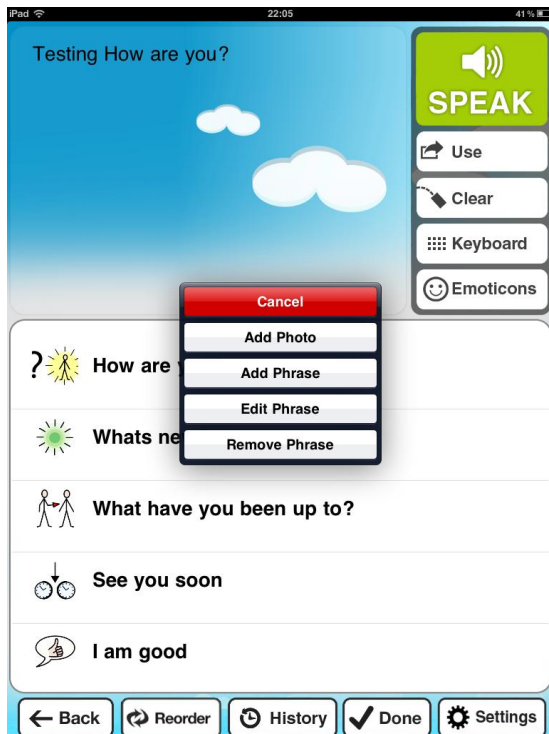


This may be useful to adjust for spelling errors or to rephrase. If you want to remove a phrase, you select “remove phrase”.

Remember after making any changes choose “Done” to come out of editing mode.

CHAPTER 12

HISTORY

History pages bring up the most recent phrases or words used.

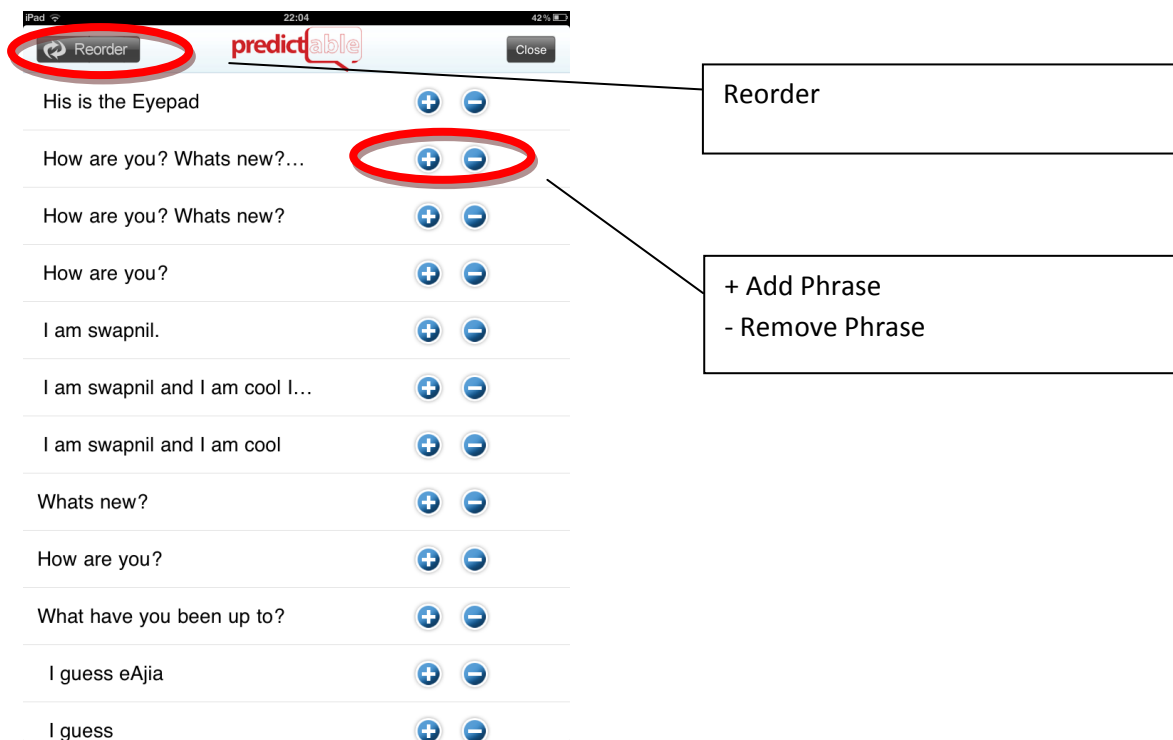

These are listed in order of most recent first. The list will display the most recent 25 items. You can select which phrase/word you want to speak or publish. This is scannable to use with your switch (or in Touch Anywhere mode). In scanning mode, the number of phrases listed is limited to that which fits onto the iDevice's screen. By selecting a phrase it will speak it directly and bring it up to the message window.

You can also choose to Add the phrase by choosing the blue + sign. You can remove the phrase by tapping the blue – sign. You can reorder the phrases by choosing the Reorder button at the top of the screen.

Note the saving of items in History can be turned off. Once turned off it will not save any further messages in History.


CHAPTER 13

BACK UPS

There are two ways you can back up – backing up directly on your iPad/iPod/iPhone; or backing up to a computer via a Wi-Fi connection.

Backing up on the iDevice:


1. Go to Settings in Predictable
2. Select Backup now

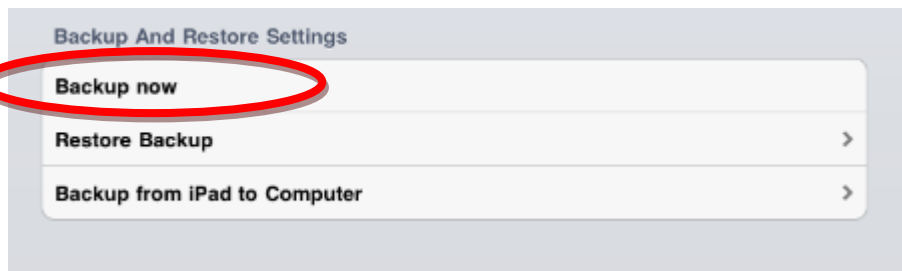

A message will pop up confirming the backup is complete.

You can then restore to a former version, as saved in backups.


1. Go to Restore Backup

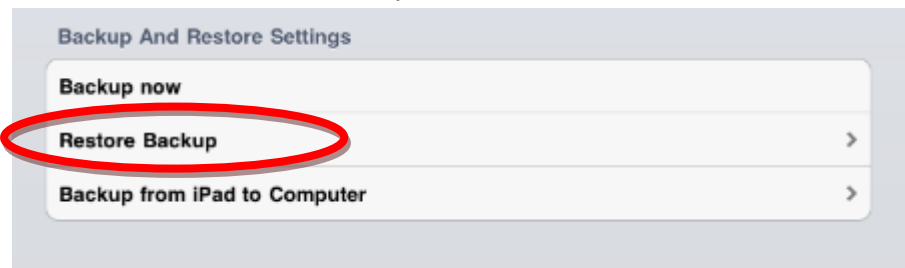

2. Choose a backup to restore to. Note they are date and time listed. You have three options:




- i) email (to either the default "to" email address you've set up or to a different email address)

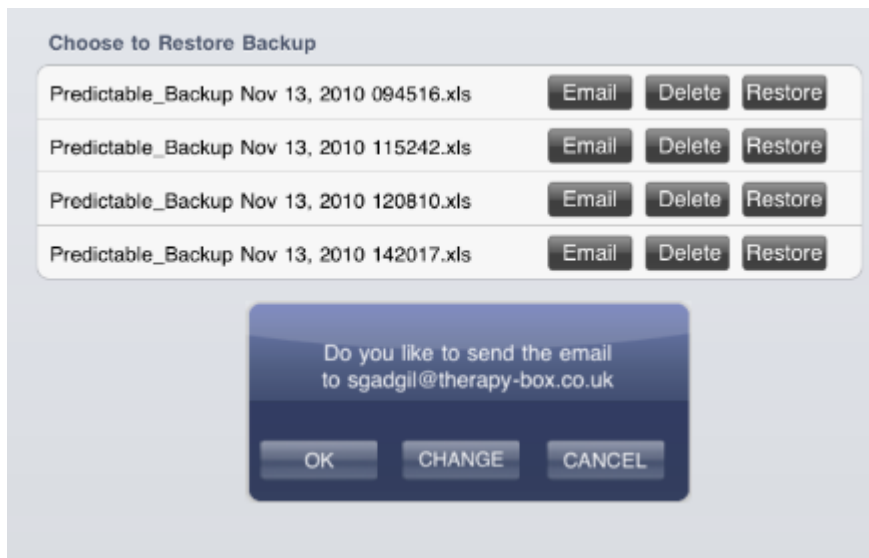


ii) delete

iii) restore


You can also backup from your iPad/iPod/iPhone to your computer. You will need Wi-Fi to do this.

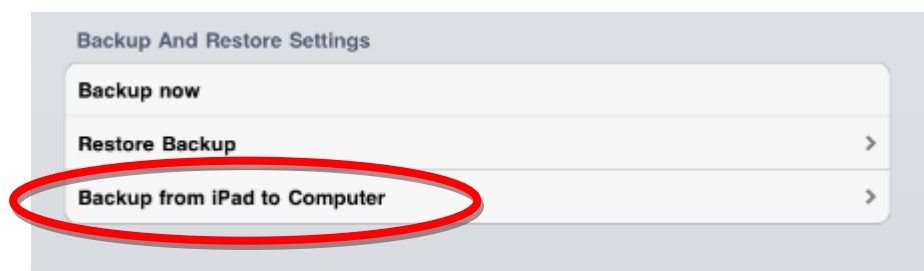

This screen will then appear, and you need to select “on” to connect to the computer. Make sure the computer and the iDevice both have Wi-Fi turned on, and are sharing the same Wi-Fi network.

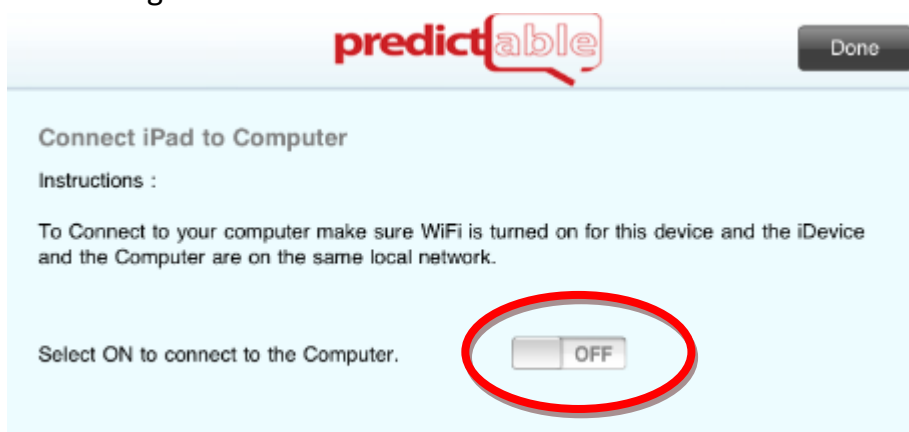

Then you will be asked to enter a web address into your browser. You can only use Safari or Firefox web browsers. You cannot use Internet Explorer or other web browsers. The web address to enter in your browser will be a string of digits.

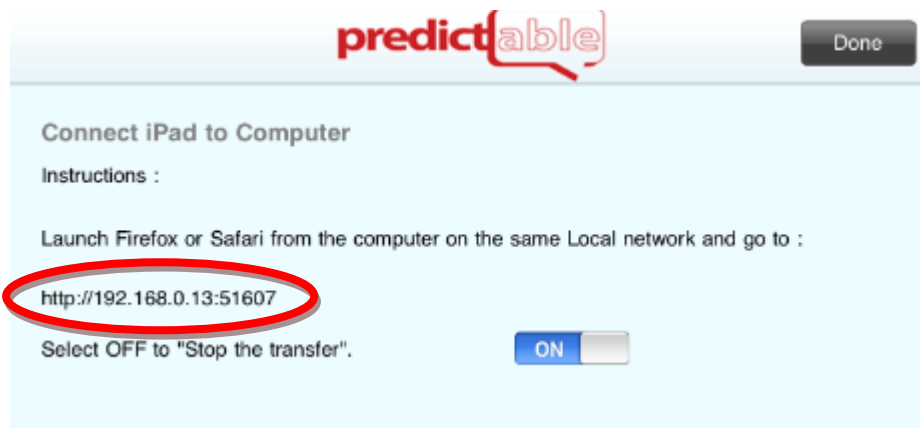

The following webpage will then come on screen:

Welcome to Predictable back up management website, your one stop shop to sending and receiving Predictable back-ups!

Useful Information

- You can view back ups in Microsoft Excel.
 - You can simply Add, Delete or Edit the phrases from any back up! Just remember to save the file before resending to Predictable on the same or any other iDevice.
 - After downloading each back up, you can rename the files for easy filling.
 - Any re-named back up simply be sent back to Predictable app on the same or other iDevice
- Let's get started...

Send a back ups from this to computer to predictable

Upload backup file from your Machine :

Download back ups from Predictable to this computer

..

[Predictable Backup Nov 13, 2010 094516.xls](#) (3.2 Kb, 2010-11-13 09:45:16 +0000)

[Predictable Backup Nov 13, 2010 115242.xls](#) (3.2 Kb, 2010-11-13 11:52:42 +0000)

This gives two options:

- 1) Sending a back up from this computer to Predictable
- 2) Download a back up from Predictable to this computer

To Send a back up from this computer to Predictable:

1. Upload a backup file already saved on your computer. You can browse for a file. When you find it, select "submit"

Welcome to Predictable back up management website, your one stop shop to sending and receiving Predictable back-ups!

Useful Information

- You can view back ups in Microsoft Excel.
 - You can simply Add, Delete or Edit the phrases from any back up! Just remember to save the file before resending to Predictable on the same or any other iDevice.
 - After downloading each back up, you can rename the files for easy filling.
 - Any re-named back up simply be sent back to Predictable app on the same or other iDevice
- Let's get started...

Send a back ups from this to computer to predictable

Upload backup file from your Machine :

Download back ups from Predictable to this computer

..

[Predictable Backup Nov 13, 2010 094516.xls](#) (3.2 Kb, 2010-11-13 09:45:16 +0000)

[Predictable Backup Nov 13, 2010 115242.xls](#) (3.2 Kb, 2010-11-13 11:52:42 +0000)

To download a back up from Predictable to the computer

1. Click on the backups which you have sent from the iDevice earlier

Welcome to Predictable back up management website, your one stop shop to sending and receiving Predictable back-ups!

Useful Information

- You can view back ups in Microsoft Excel.
 - You can simply Add, Delete or Edit the phrases from any back up! Just remember to save the file before resending to Predictable on the same or any other iDevice.
 - After downloading each back up, you can rename the files for easy filling.
 - Any re-named back up simply be sent back to Predictable app on the same or other iDevice
- Let's get started...

Send a back ups from this to computer to predictable

Upload backup file from your Machine :

Download back ups from Predictable to this computer

..

[Predictable Backup Nov 13, 2010 094516.xls](#) (3.2 Kb, 2010-11-13 09:45:16 +0000)

[Predictable Backup Nov 13, 2010 115242.xls](#) (3.2 Kb, 2010-11-13 11:52:42 +0000)

2. You will be presented with a page which gives you information on settings and the phrases stored in My Phrases. It looks like this. Please note "rowend" indicates an empty category.


Predictable							
Settings							
Voice Selected	Alan	Change Option are: Alan, Allison, Elizabeth, Grace, Kate, Simon, Steven, Susan					ROWEND
Pitch Rate Selected	Normal	Change Option are: Deepest, Deeper, Deep, Normal, High, Higher, Highest					ROWEND
Speech Rate Selected	Normal	Change Option are: Very Slow, Slow, Normal, Fast, Very Fast					ROWEND
Phrases							
CHAT	How are you?	Whats new?	What have you been up to?	See you soon	I am good	Good to see you	ROWEND
HELP	I need help	I am hurting	I am thirsty	I need to go to the toilet	I want to go to bed	ROWEND	
FOOD & DRINK	Coffee with milk please	A cup of tea please	ROWEND				
QUESTIONS	ROWEND						
AT HOME	ROWEND						
OUT & ABOUT	ROWEND						
THINGS I LIKE	ROWEND						
MY FAMILY	ROWEND						
MY FRIENDS	ROWEND						
STUFF I SAY ALL THE TIME	ROWEND						
FEELINGS	ROWEND						
MY STUFF	ROWEND						

CHAPTER 14

EMAIL

You can send an email from Predictable. Remember you will need a WiFi or 3G connection to be able to do this.


After typing in your message window the content of the email, you will choose “Use” which is in the keys on the right of the message window. Then choose Email.

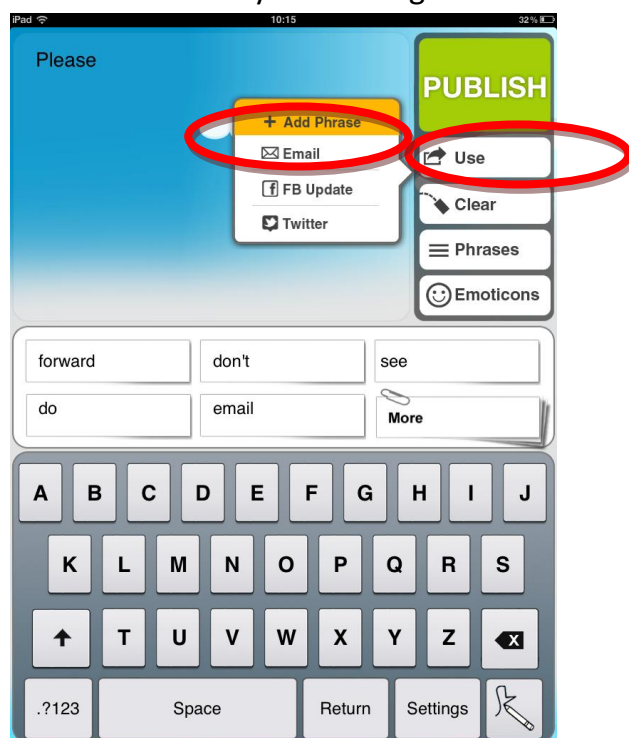

This is accessible via direct touch or scanning.

You will be then asked if you would like to send the email to the address entered in the settings page, or if you would like to change. To change, you will need to either

enter the email address directly or choose from your email contact list imported via iTunes..

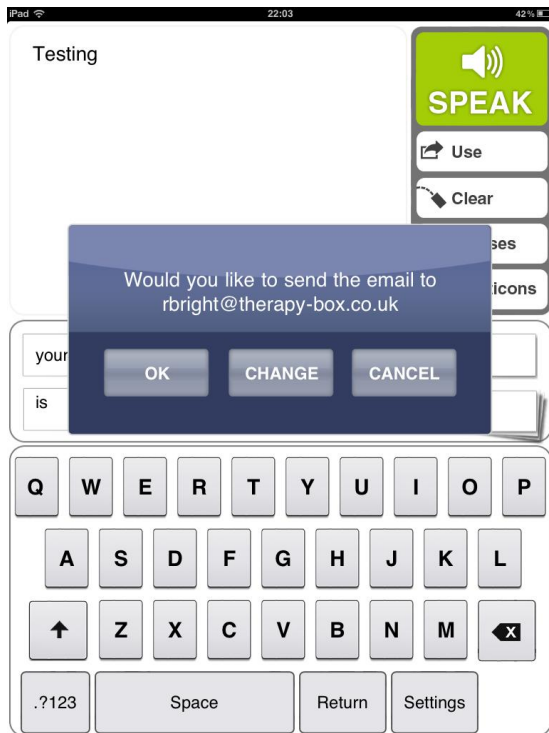

After clicking “ok” the email will be sent.

To set up the email system, you need to:

1. Go to Settings on your iPad/iPhone/iPad Touch



2. Scroll down to Predictable



Go to Mail Settings. Enter a “from” email address with the password that works for that address. Choose from the following email hosts:

Yahoo

Hotmail

Google

Mobile Me

Others

You will then need to enter your email address and password to be saved in the settings. If you have an email host other than the four listed above, you would need to enter the host smtp address.


You can also enter a default To email address if there is one email address that the person using Predictable would like to primarily send emails to. Within the app, they can choose this default address or search through their other contacts.

CHAPTER 15

FACEBOOK AND TWITTER

Please note that updating your Facebook status or tweeting is only possible in direct access mode. It is not currently available in switch or scanning modes. You will be

asked to login with your Facebook or Twitter user name and password.

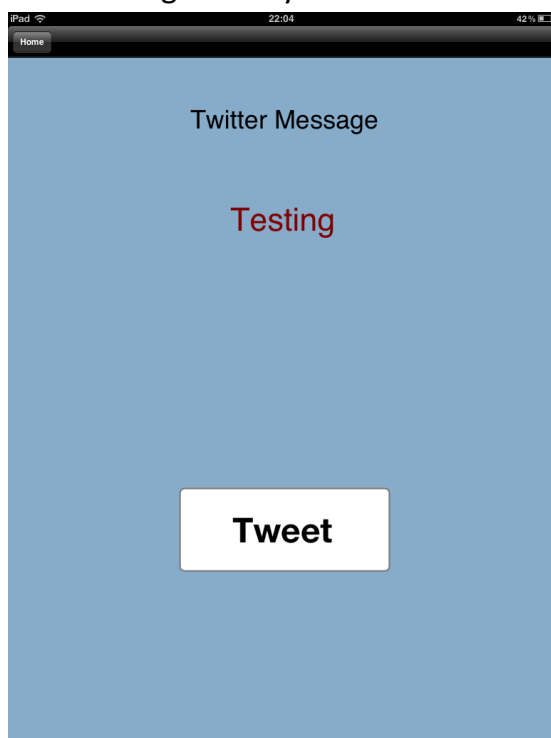

CHAPTER 16

SENDING AN SMS

SMS can only be sent from the iPhone. From the settings menu at the bottom of the page, select the SMS option. Any message you had just created (and not deleted) in the message window will then transfer to the content of your SMS. Now you can enter the name from your Contacts on your iPhone or enter the telephone number and select Send.

Please note that while it is possible to create the original message via scanning or switch mode, it is not possible to send the SMS via scanning or switch.


CHAPTER 17

EMOTICONS

There is an option for you to add Emoticons to your spoken message. Emoticons are sounds including cough, cry, eh, kiss, mmm, oh, laugh, throat clear, whistle and yawn. These can be added at the end or within a spoken message or used in isolation. If the message window is clear and you choose an emoticon it will play this sound on its own.

If used when compiling a message the emoticon will appear in the top right corner of the message window.

The emoticon will be saved in History and Category Phrases if used.

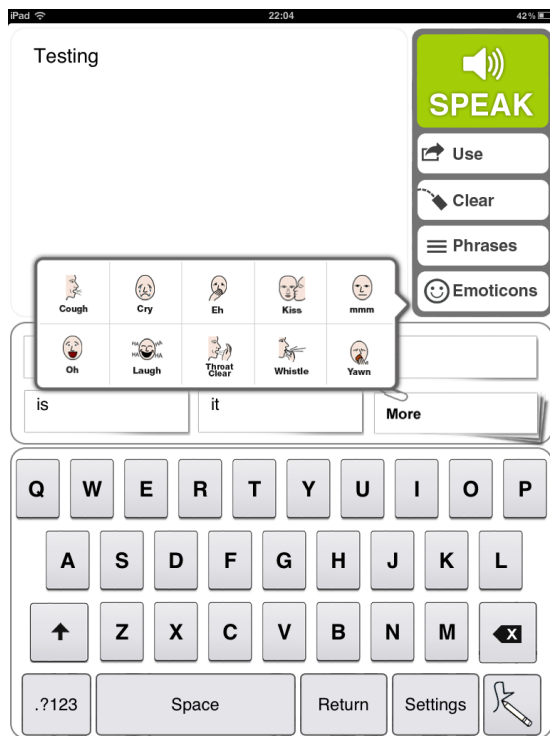


This feature can also be turned off by going to the Settings on the iDevice (under Predictable settings).

Note some of the Loquendo voices do not have all of the Emoticons available for them.

CHAPTER 18


USING HAND WRITING INPUT

The option of using handwriting input is now possible for direct access only. You can write a letter or word and then select it to go in the message window.

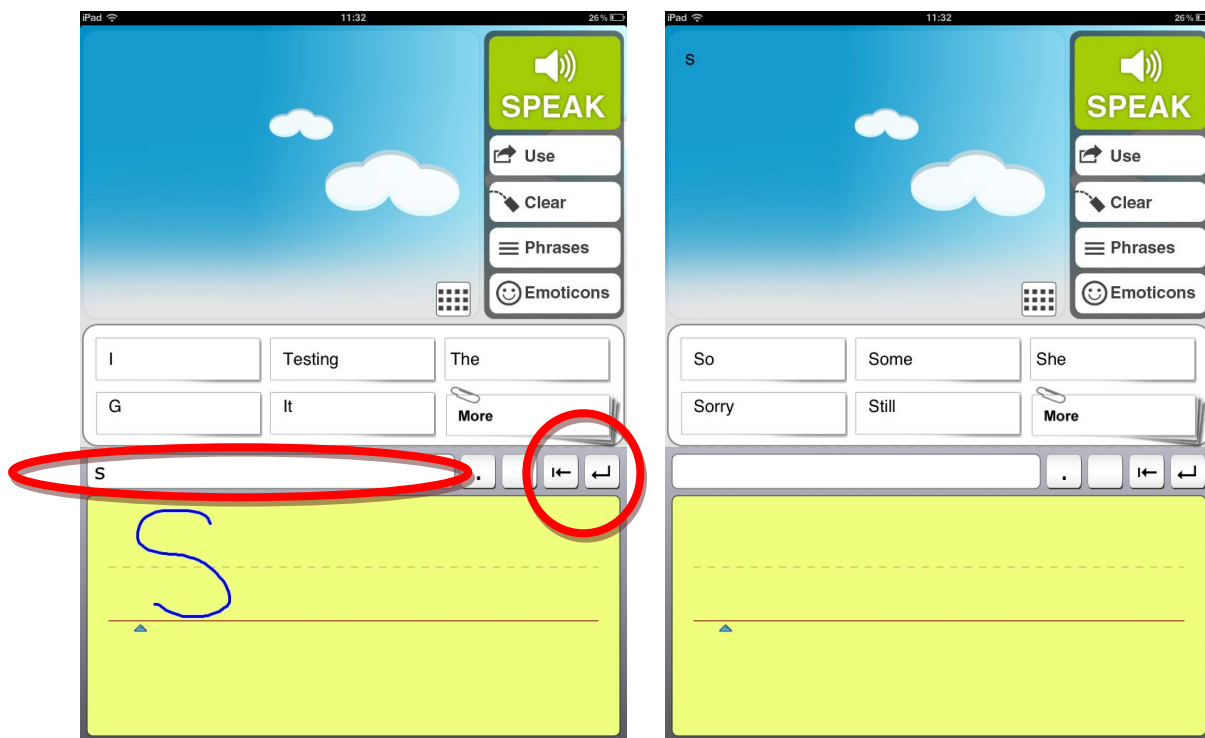

It acts like any other keyboard entry so you can use it with the prediction options too.

First write the letter/s. Then check that the note pad is accurately recognising it by looking at the small window directly above the note pad.

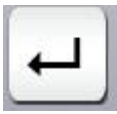If it is correct choose the return arrow key  to select it. To erase the writing

choose the back arrow key . To return to the keyboard, select the keyboard

key on the message window space.



CHAPTER 19

FEATURE ACCESSIBILITY – VIA DEVICE AND ACCESS METHOD

Features	iPod Touch	iPhone	iPad
SMS	No	Yes in direct access only	No

Facebook Status Update	Yes with Wifi in direct access only	Yes with Wifi or 3G in direct access only	Yes with Wifi in direct access only
Email	Yes with Wifi in all access modes	Yes with Wifi or 3G in all access modes	Yes with Wifi in all access modes
Twitter	Yes with Wifi in direct access only	Yes with Wifi or 3G in direct access only	Yes with Wifi in direct access only.
Note Pad	Yes in direct access only	Yes in direct access only	Yes in direct access only